



## **STERIGERMS 12 Litres User's Manual**



**STERIGERMS International LLP**

**Unit 3 – 8th floor Ellerman LLP**

**12 – 20 Camomile street**

**LONDON EC3A 7PT**

**Phone: +33 (0)4.90.11.37.35 - Fax: +33 (0)4.90.11.37.30**

**VAT number : GB 8722853 06**




## TABLE OF CONTENTS

<b>1. OVERVIEW .....</b>	<b>7</b>
<b>1.1. Introduction (Figure 1) .....</b>	<b>7</b>
<b>1.2. Characteristics.....</b>	<b>8</b>
1.2.1. Physical characteristics .....	8
1.2.2. Electrical characteristics.....	8
<b>2. OPERATION .....</b>	<b>9</b>
<b>2.1. Connection.....</b>	<b>9</b>
<b>2.2. Power-up .....</b>	<b>9</b>
2.2.1. Power-up after shutdown .....	9
2.2.2. Standard power-up.....	10
<b>2.3. Exploitation cycle.....</b>	<b>12</b>
<b>2.4. Case of wrong password.....</b>	<b>21</b>
<b>3. SHUTDOWN .....</b>	<b>21</b>
<b>3.1. Standard shutdown.....</b>	<b>21</b>
<b>3.2. Complete shutdown .....</b>	<b>21</b>
<b>3.3. Emergency shutdown .....</b>	<b>22</b>
<b>4. ALARMS .....</b>	<b>23</b>
<b>5. MAINTENANCE .....</b>	<b>25</b>
<b>5.1. Table of periodic operations .....</b>	<b>25</b>
<b>5.2. Cleaning of the treatment compartment .....</b>	<b>25</b>
<b>5.3. Control of oil level .....</b>	<b>26</b>
<b>5.4. Control of cooling liquid level.....</b>	<b>27</b>
<b>5.5. Exchange of the cover gasket.....</b>	<b>28</b>
<b>5.6. Exchange of the particle filter .....</b>	<b>28</b>
<b>6. WARNING AND SAFETY SIGNS.....</b>	<b>31</b>
<b>7. HUMAN – MACHINE DIALOG .....</b>	<b>33</b>
<b>8. WARRANTY.....</b>	<b>53</b>



## LIST OF HUMAN – MACHINE INTERFACE SCREENS

TABLE OF CONTENTS.....	3
Window 1 - Welcome Please press the “Enter” button .....	33
Window 2 - Select a user And type your password .....	34
Window 3 – Correct password .....	35
Window 4 – Wrong password Please wait and type it again.....	36
Window 5 - Warning (first window).....	37
Window 6 – Warning (second window) .....	38
Window 7 - Warning (third window).....	39
Window 8 – Select an operation.....	40
Window 9 - Descent of the piston .....	41
Window 10 - Ejection of the compressed waste product.....	42
Window 11 - Closing Do you really want to shut down the device .....	43
Window 12 – Press Start to Launch the cycle .....	44
Window 13 - Press Confirm to launch the cycle .....	45
Window 14 – Start compression .....	46
Window 15 - Heating.....	47
Window 16 – Start threshold .....	48
Window 17 - Cooling .....	49
Window 18 – End of cycle.....	50
Window 19 – Cleaning of the machine .....	51
Window 20 - Machine error Only administrator can enter.....	52

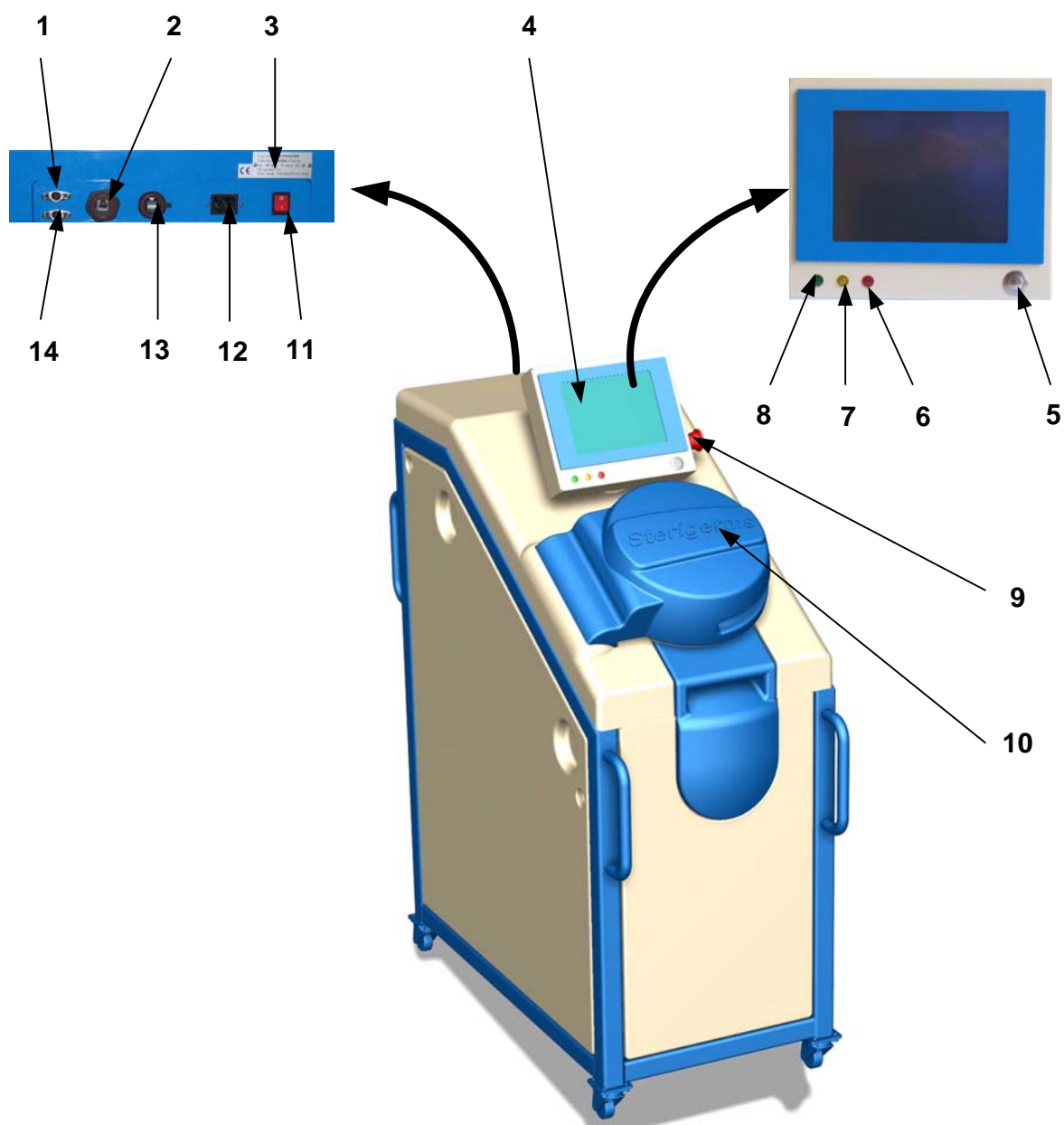

**Figure 1 - STERIGERMS**

# 1. OVERVIEW

## 1.1. Introduction (Figure 1)

Thank you for choosing STERIGERMS for the decontamination of your IMW (*Infectious Medical Waste*).

STERIGERMS ensures treatment of Infectious Medical Waste, such as dressing, yellow containers used for disposal of needles and of stinging and cutting objects, syringes, diapers, probes, blood bags, catheters, bandages, gloves ...

Combining volume reduction and decontamination by heating, STERIGERMS allows turning potentially infectious waste into a flat-shaped compressed waste product disposable as normal household waste.

STERIGERMS is easy to use thanks to its touch screen and easy to operate, as it only requires a 16-Ampere power supply. STERIGERMS can be operated close to the location of origin of the waste, thus limiting the risks of contamination implied by the transport and handling of yellow bags.

STERIGERMS is composed of the following elements:

- A locking cover (10) allowing access to the waste treatment compartment
- A control console (4) composed of:
  - a touch screen allowing human-machine dialog
  - a red indicator light (6): when lit, signals a major incident
  - a yellow indicator light (7): when lit, signals a minor incident without any impacts on the on-going cycle
  - a green indicator light (8): when lit, indicates correct operation
  - a pushbutton (5) to power up the control console
  - an emergency knock-off switch (9) to switch off mains power supply, equipped with a removable key to unlock the cover
- At the rear, a connection panel with the following connections:
  - a nameplate (3)
  - a main power-up switch (11) with light indication
  - a socket (12) to connect an external power supply
  - four connection elements not used for standard operation:
    - a USB connector (13)
    - a socket (14) to connect a keyboard
    - a socket (1) to connect a mouse
    - a socket (2) for an Ethernet link.

**1.2. Characteristics**

**1.2.1. Physical characteristics**

LENGTH (mm)	DEPTH (mm)	HEIGHT (mm)	WEIGHT (kg)
460	1 050	1 120	140

**1.2.2. Electrical characteristics**

- Power supply tension : 230V
- Frequency : 50/60Hz



## 2. OPERATION

**Caution:** Any identified malfunction of STERIGERMS is displayed on the banner of each window and must be reported to the administrator.

### 2.1. Connection

**Caution:** Before connecting a device, always check that the power supply is of type “16A - 230V” with earth connection (alternative current, 16 Amperes).

Connect the power cable to the socket (1, Figure 2) of STERIGERMS and to the external power supply.

### 2.2. Power-up


**Caution:** When powering up the equipment for the first time, contact the administrator to unlock STERIGERMS.

There are two ways of powering up STERIGERMS:

- Power-up after shutdown
- Standard power-up.

#### 2.2.1. Power-up after shutdown

- You must perform this type of power-up after a complete shutdown of the device due to a technical issue or a transfer of the device to another location.

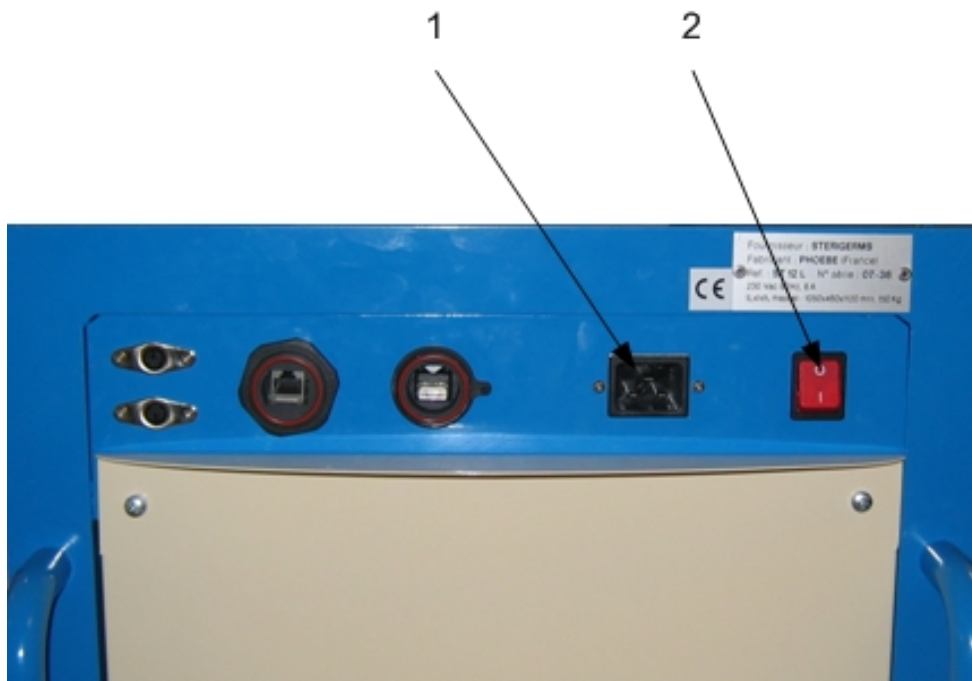


**Figure 2 – Rear side of STERIGERMS**

- Check that STERIGERMS is connected to the power supply network (§2.1).
- Ensure that the emergency knock-off switch (9, Figure 1) located at the right-hand side of the STERIGERMS control console is directed outwards; if not, bring it back to its original position using the key and then retrieve the key.

- Set the main power-up switch (2, Figure 2) located at the rear of STERIGERMS to the “ON” position and check that the integrated light indicator comes on.
- Push the control console pushbutton (4, Figure 3).

### 2.2.2. Standard power-up

- Ensure that the emergency knock-off switch (9, Figure 1) located at the right-hand side of the STERIGERMS control console is directed outwards ; if not, put it back into its original position using the key and then retrieve the key.
- Check that the main power-up switch (2, Figure 2), located at the rear of STERIGERMS is in the “ON” position and that the integrated indicator light is on.
- Push the control console pushbutton (4, Figure 3).

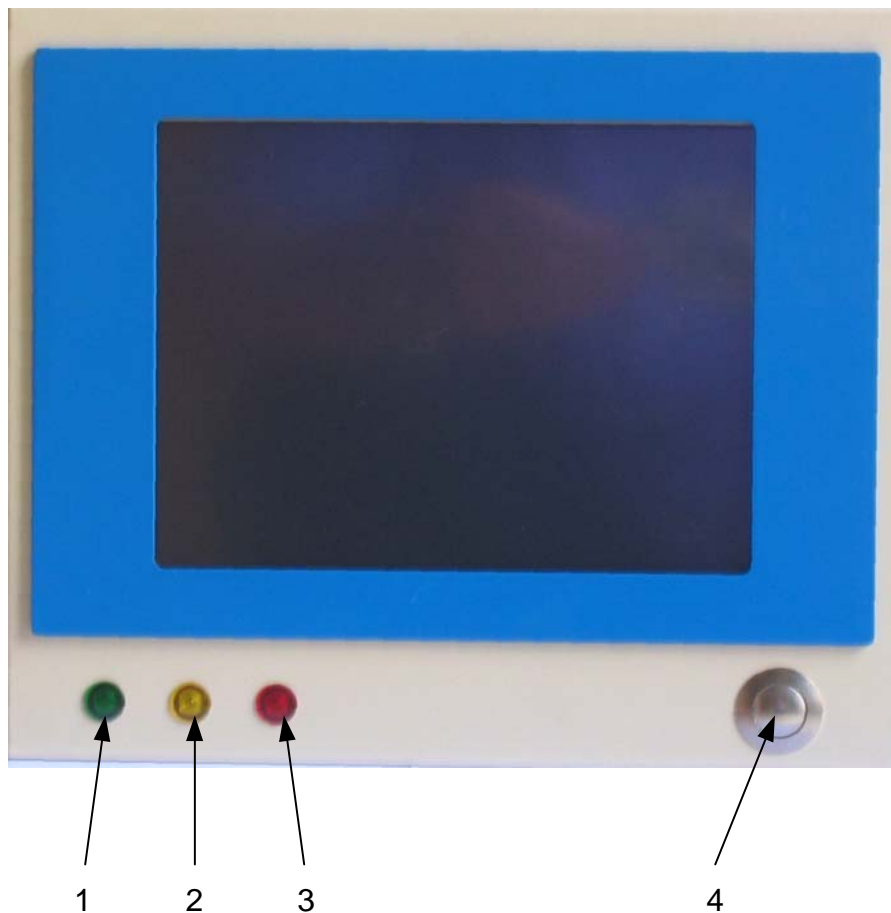

**Figure 3 – Control console**

The control console red (3, Figure 3) and yellow (2, Figure 3) indicator lights come on for about thirty seconds (time required for initialization of Windows and Sterisoft software) and go out. Then, the green indicator light (1, Figure 3) comes on and the Welcome page appears.



Caution: In case of problem (if the red or yellow indicator light(s) stay(s) on), stop the power-up procedure and contact the administrator. If the yellow indicator light comes on during a cycle, the cycle can be completed anyway.

## 2.3. Exploitation cycle

- On the screen, select the “Enter” button. The **‘Select a user And type your password’** window appears:

The screenshot shows a software window titled "StériSoft 12L V4.031c". The main heading inside the window is "Select a user And type your password". The interface is divided into several sections:

- User Selection:** A vertical column of five buttons labeled "User 1", "User 2", "User 3", "User 4", and "Administrator".
- Keypad:** A 3x3 grid of buttons. The first two rows contain digits 1-3, 4-6, and 7-9. The third row contains "Del.", "0", and "Reset".
- Input and Action:** A text input field at the top right, a "Validate" button below it, and an "Exit" button at the bottom right.
- Navigation:** A "Quit" button at the bottom left and an "Enter" button at the bottom center.
- Status Bar:** A box at the bottom right displaying the date and time: "23/05/2007 17:40:38".

- Select the required user. This highlights the corresponding button.
- Using the keyboard, type in the corresponding 4-digit password and press the “Validate” button. The **‘Correct password’** window appears. If the typed password is incorrect, the **‘Wrong password’** window is displayed. In that case, refer to § 2.4.

StériSoft 12L V4.031c

## Correct password

User 1	1	2	3	XXXX
User 2	4	5	6	Validate
User 3	7	8	9	<b>Correct password</b>
User 4	Del.	0	Reset	Exit
Administrator				

Quit
Enter

23/05/2007  
17:44:36

- On the screen, select the “Enter” button. The first ‘**Warning**’ window appears.

StériSoft 12L V4.031c

## Warning


The use of the STERIBAG bag with the STERIGERMS device is compulsory. The STERIBAG bag have imperatively to be used exclusively from any other type of content, considering that they are the only ones warranting the sterilized results targeted, provided all the other conditions of process are implemented. Consequently the seller wave all responsibilities with regards to any damages of any kind as a result of non-respect of this particular provision such being the basis of the Agreement delivered by the relevant Authority. Further the seller warranties to have all buyers or whole-sellers within the chain of distribution implement this provision.

Validate


Quit
Exit

23/05/2007  
17:43:53

- On the screen, select the “Validate” button. The ‘**Select an operation**’ window appears.

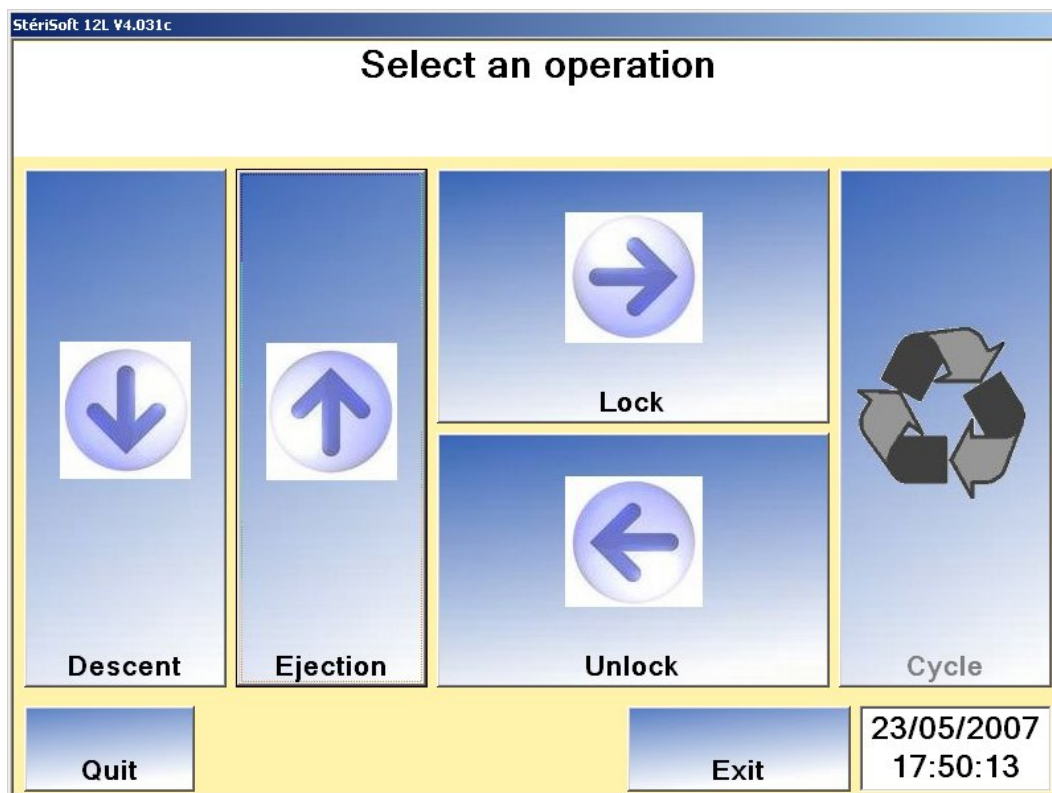

- On the screen, select the “Descent” button to move down the piston. The following window appears:

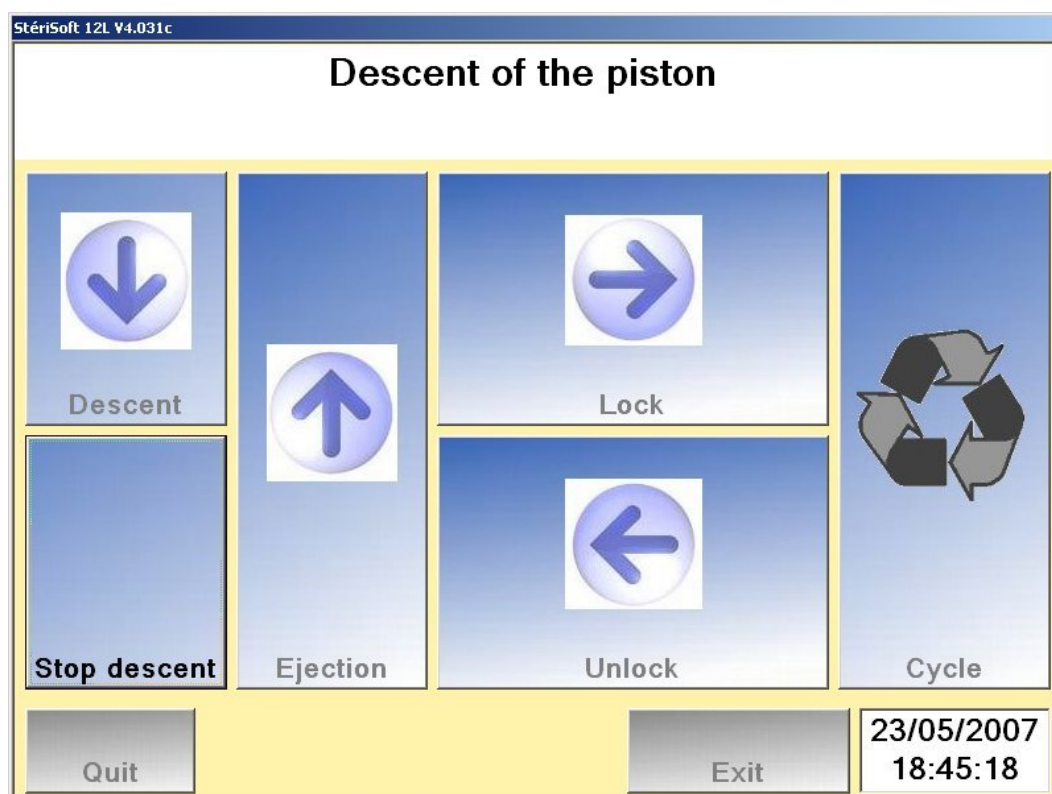

- Insert the waste bag into the Steribag bag for treatment,
- Insert the Steribag bag into the STERIGERMS at the bottom of the waste treatment compartment so that it does not overflow the edges of the container.




- 1) Fold the lower part of the steribag
- 2) Fold the top part of the steribag




- 3) Fold down the two tips thus formed along the barrel

**Figure 4 – Waste bag closing procedure**


- Close the cover of STERIGERMS.
- On the screen, select the “Lock” button (locking lasts about 5 seconds).
- On the screen, select the “Cycle” button, the second ‘**Warning**’ window appears.

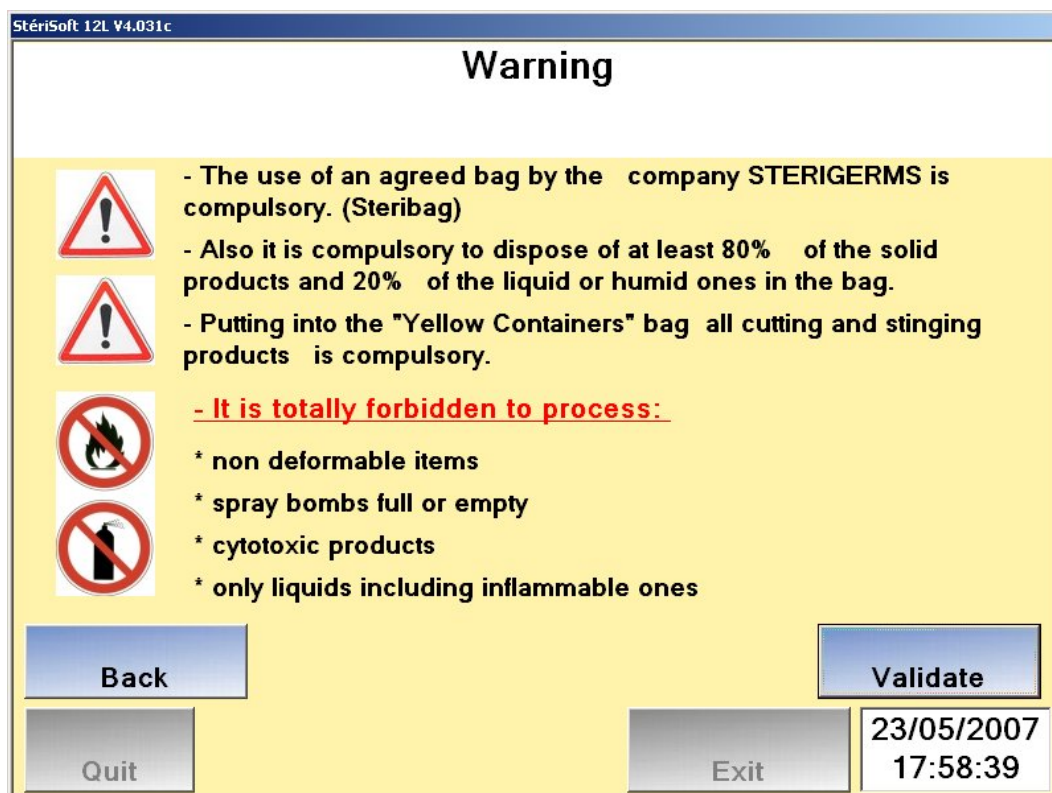

- On the screen, select the “Validate” button, the third ‘**Warning**’ window appears.

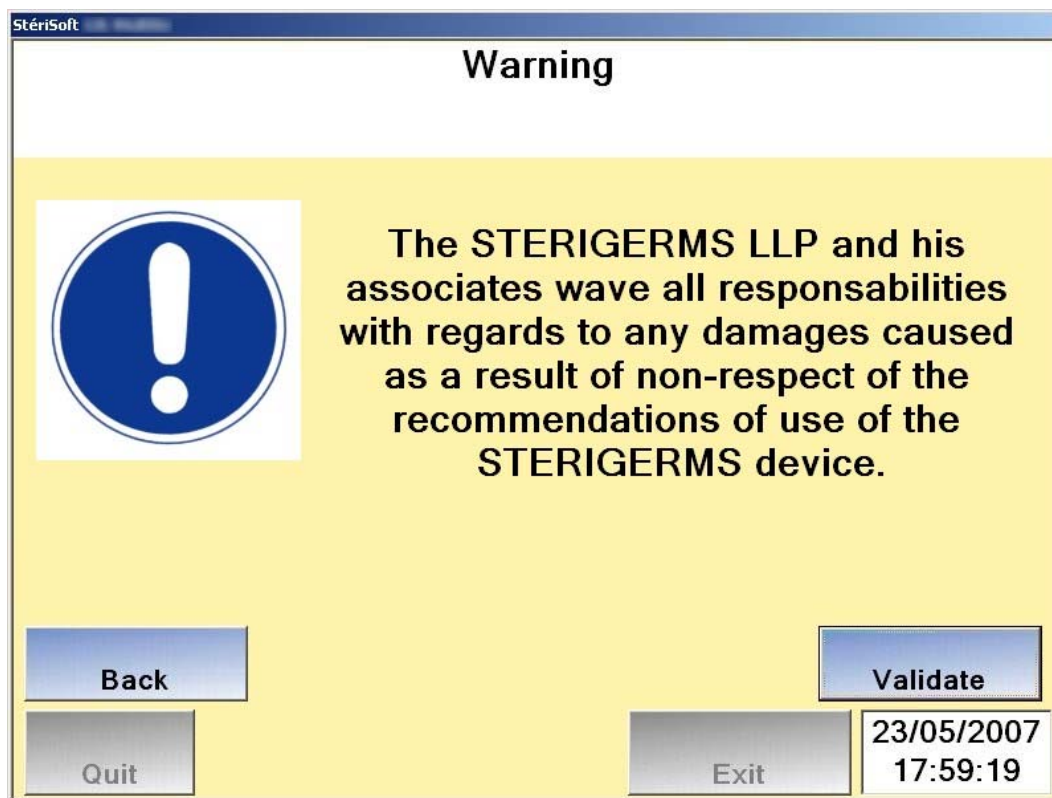

- On the screen, select the “Validate” button to access the ‘**Press Start to launch the cycle**’ window.

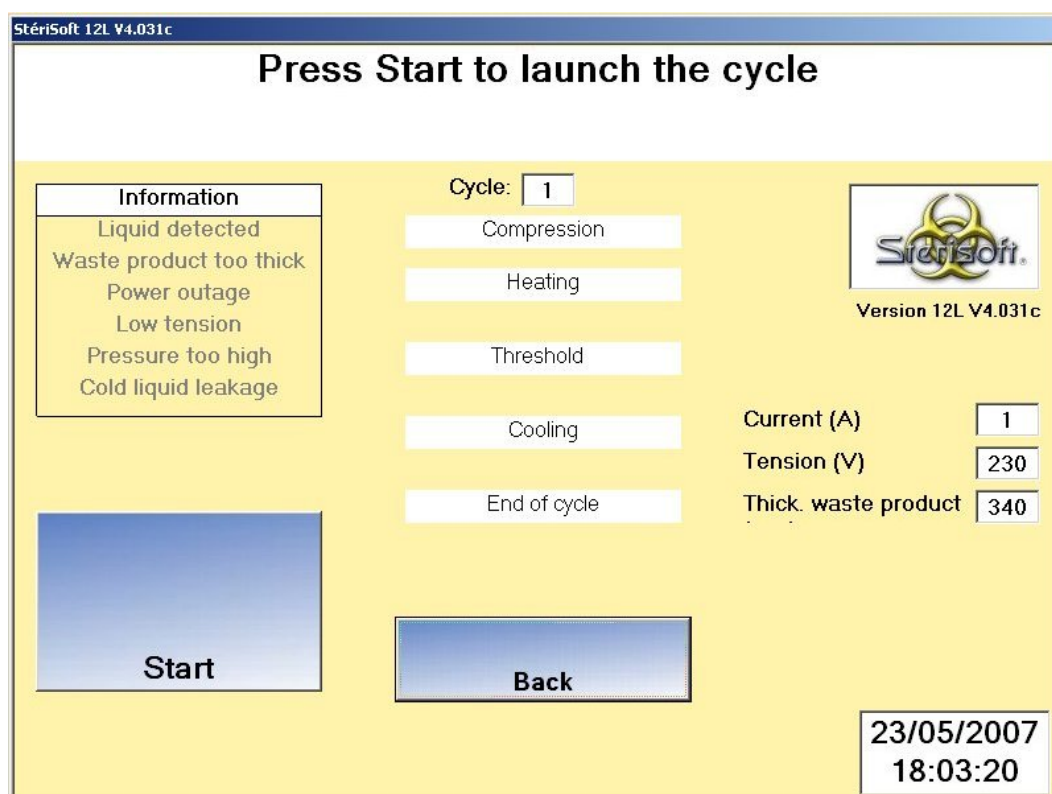

- On the screen, select the “Start” button to access the ‘**Press Confirm to launch the cycle**’ window.



StériSoft 12L V4.031c

## Press Confirm to launch the cycle

Information
Liquid detected
Waste product too thick
Power outage
Low tension
Pressure too high
Cold liquid leakage

Cycle:


Compression

Heating


Threshold

Cooling

End of cycle



Version 12L V4.031c



**Confirm**

**Back**

Current (A)

Tension (V)

Thick. waste product

**23/05/2007**  
**18:04:48**

- On the screen, select the “Confirm” button to launch the cycle.

During the cycle, the following windows automatically appear to visualize the various stages of the ongoing cycle (compression, heating, treatment, cooling, end of cycle).

StériSoft 12L V4.031c

## Starting compression

Information
Liquid detected
Waste product too thick
Power outage
Low tension
Pressure too high
Cold liquid leakage

Cycle:


Compression

Heating


Threshold

Cooling

End of cycle



Version 12L V4.031c



Current (A)

Tension (V)

Thick. waste product

**23/05/2007**  
**18:30:21**

StériSoft 12L V4.031c

## Heating

Information

Liquid detected  
Waste product too thick  
Power outage  
Low tension  
Pressure too high  
Cold liquid leakage

Cycle: 1


Compression

Heating  
23/05/2007 18:33:34


Threshold

Cooling

End of cycle


  
Version 12L V4.031c

Current (A) 1  
Tension (V) 230  
Thick. waste product 53



23/05/2007  
18:33:37

StériSoft 12L V4.031c

## Threshold reached

Information

Liquid detected  
Waste product too thick  
Power outage  
Low tension  
Pressure too high  
Cold liquid leakage

Cycle: 1


Compression

Heating  
23/05/2007 18:35:55


Threshold remain 12  
23/05/2007 18:35:55

Cooling

End of cycle


  
Version 12L V4.031c

Current (A) 1  
Tension (V) 230  
Thick. waste product 53



23/05/2007  
18:36:01

Stérisoft 12L V4.031c

## Cooling

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


Compression

Heating  
23/05/2007 18:37:30

Threshold  
23/05/2007 18:37:30

Cooling  
23/05/2007 18:37:30

End of cycle




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product



23/05/2007  
18:37:37

Stérisoft 12L V4.031c

## End of Cycle

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


Compression

Heating  
23/05/2007 18:41:03

Threshold  
23/05/2007 18:41:03

Cooling  
23/05/2007 18:41:03

End of cycle  
23/05/2007 18:41:03




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product




Waste has been processed

Opening device


23/05/2007  
18:41:08

When the cycle has completed, the message “Waste has been processed / Waste has not been processed” is displayed in the lower part of the screen.

- On the screen, select the “Open device” button to return to the ‘**Select an operation**’ window.

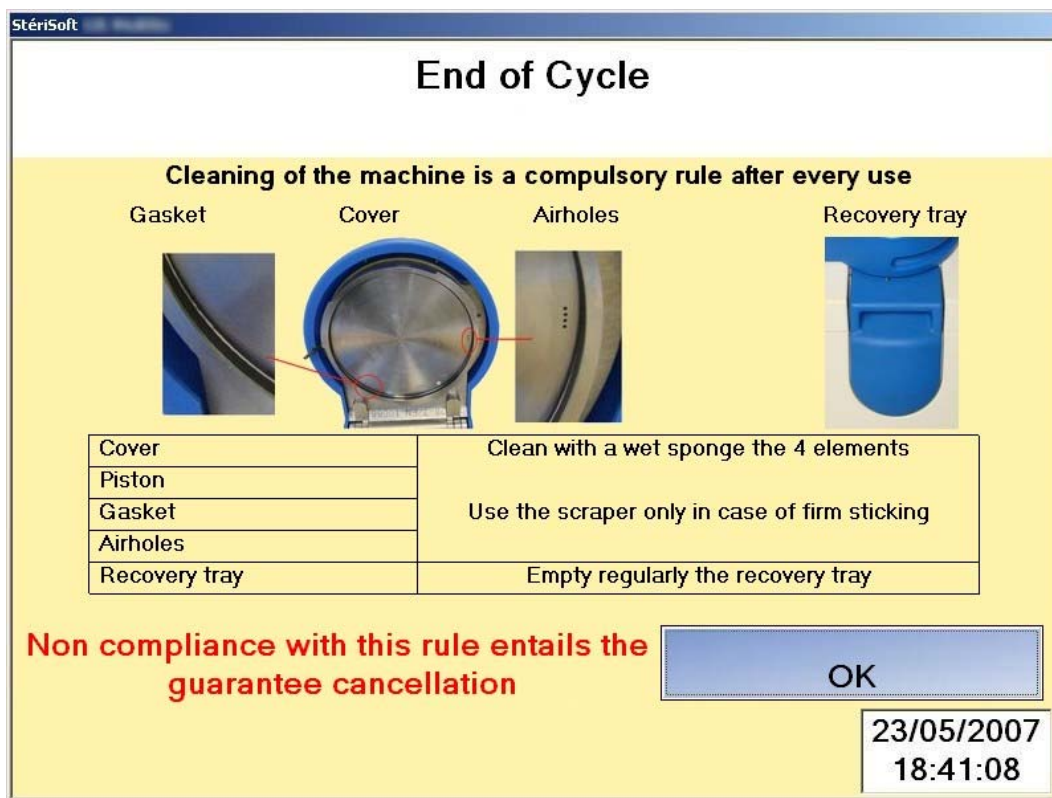

- On the screen, select the “OK” button. The “Select an operation” window appears

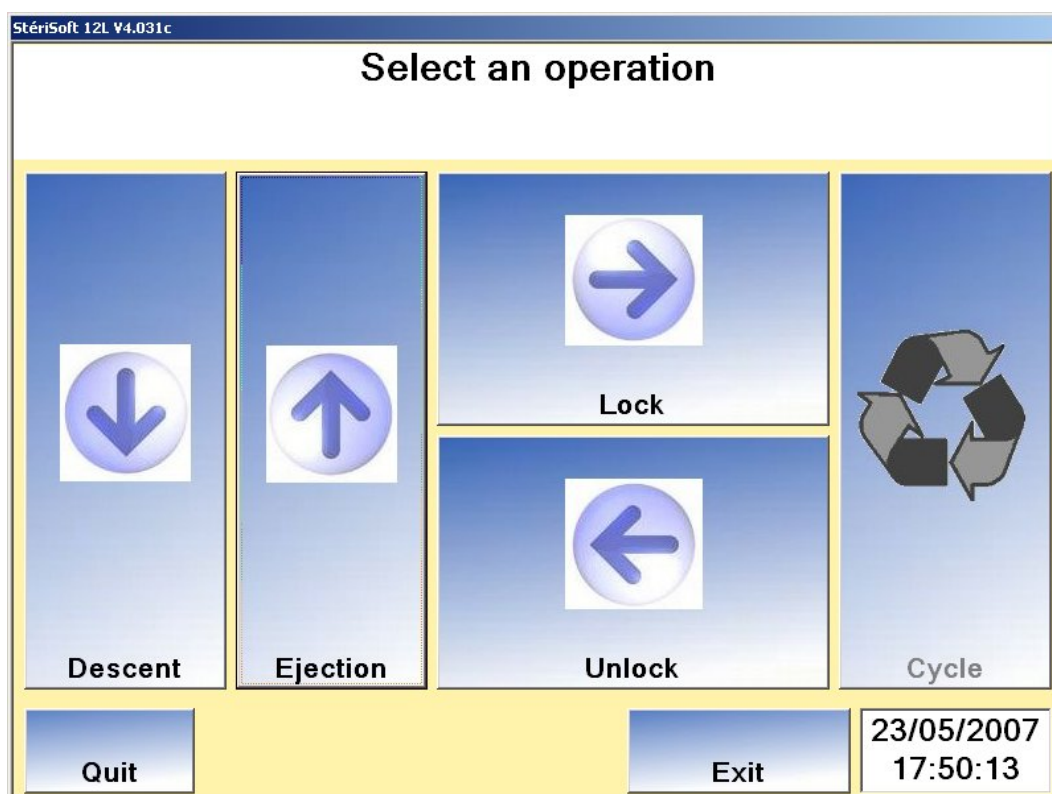


- On the screen, select the “Unlock” button.
- Open the cover of STERIGERMS.
- On the screen, select the “Eject” button to monitor the ejection of the piston.

- Retrieve the compressed waste product and discard it. You can now dispose of it as of standard household waste.

## 2.4. Case of wrong password

If the password typed is wrong, the '**Wrong password Please wait and type it again**' window appears during 4 seconds, then the '**Select a user And type your password**' window is displayed again. Perform again the procedure described in § 2.3.

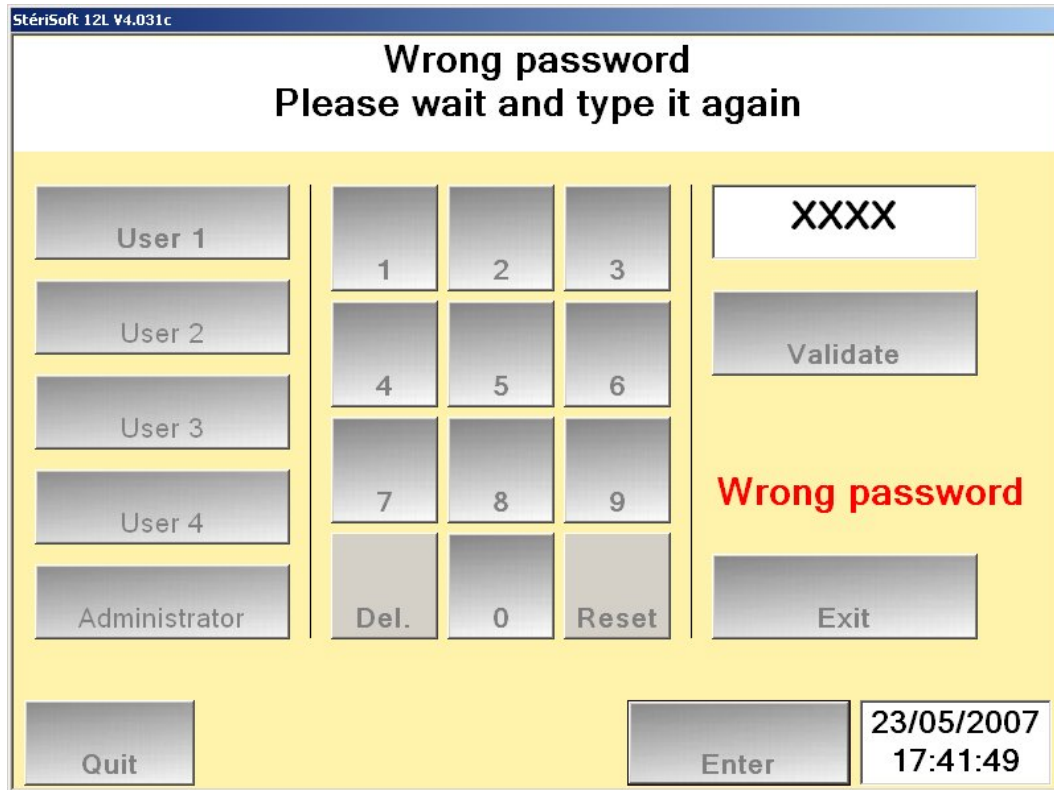

## 3. SHUTDOWN

Three types of shutdown can be performed:

- Standard shutdown after use
- Complete shutdown due to technical issue or transfer of the device to another location
- Emergency shutdown in case of serious incident

### 3.1. Standard shutdown

- On the displayed screen, select the "Quit" button and then the "YES" button of the confirmation window. Only the control console is shut down, STERIGERMS is in stand-by mode.

### 3.2. Complete shutdown

- On the displayed screen, select the "Quit" button and then the "YES" button of the confirmation window. Only the control console is shut down, STERIGERMS is in stand-by mode.
- Set the main power-up switch (2, Figure 2) located at the rear of STERIGERMS to the "Off" position. The integrated light indicator goes out.

### **3.3. Emergency shutdown**

- Push the emergency knock-off switch (9, Figure 1) located at the right-hand side of the STERIGERMS control console. This action will switch off main power supply.



## 4. ALARMS

DISPLAYED ALARMS	EXPLANATION	ACTIONS
Liquid detected	Presence of liquid under the cover	None
Waste product too thick	Thickness of compressed waste product exceeds 12 cm	None
Power outage	Once power supply has been restored, cycle is restarted from the point of interruption caused by power outage	None
Low tension	Tension of external power supply is less than 215 V	Contact the administrator
Pressure too high	Cylinder pressure too high	Contact the administrator
Lack of cooling liquid	Cooling liquid level too low	Fill up with cooling liquid
Upper banner of the screen is red (see following window)	Internal malfunction	Contact the administrator







## 5. MAINTENANCE

### 5.1. Table of periodic operations

OPÉRATIONS	FREQUENCY	PROCÉDURES
Cleaning of the treatment compartment	After each cycle	§ 5.2
Control of oil level	Twice a month	§ 5.3
Control of cooling liquid level	In case of alarm	§ 5.4
Exchange of the cover gasket	Twice a year or after production of 300 compressed waste products	§ 5.5
Exchange of the particle filter	Twice a year or after production of 300 compressed waste products	§ 5.6

### 5.2. Cleaning of the treatment compartment

- Do not use any abrasive, scrubbing or metallic sponge or sharp or cutting object for cleaning
- Clean the treatment compartment, the cover and the cover gasket surface using kitchen towel
- You can optionally use the provided scraper (Figure 5) to remove plastic bag residues on parts made of stainless steel (cover and piston)



**Figure 5 - Scraper**

- Lubricate the gasket with a standard oil that supports temperatures exceeding 145 °C
- Remove residual stains (blood, gelose...) with a sponge and regular liquid soap.

### 5.3. Control of oil level


During this operation, the piston must be kept in low position:

- Unscrew the two ¼" fastening screws (Pozidriv no 1 screwdriver) on the left flank. Tilt the flank down on your side and pull it upwards to remove it.
- Check that oil level is within the limits indicated on the oil tank.
- If oil level is correct, engage the right flank on the locating pins located at the bottom of the tank, press the flank against the chassis and tighten it using the two ¼" fastening screws (Pozidriv no 1 screwdriver).
- If oil level is too low, top up with oil (32 hydraulic synthetic oil):
  - Unscrew the two ¼" fastening screws (Pozidriv no 1 screwdriver) on the right flank. Tilt the flank down on your side and pull it upwards to remove it.
  - Remove the oil tank cap (1).
  - Insert a flexible funnel of the same type as the one shown on Figure 6.



**Figure 6 – Flexible funnel**

- Top up with oil (with a level slightly below the maximum indicated on the tank).
- Remove the flexible funnel.
- Replace the oil tank cap (1).

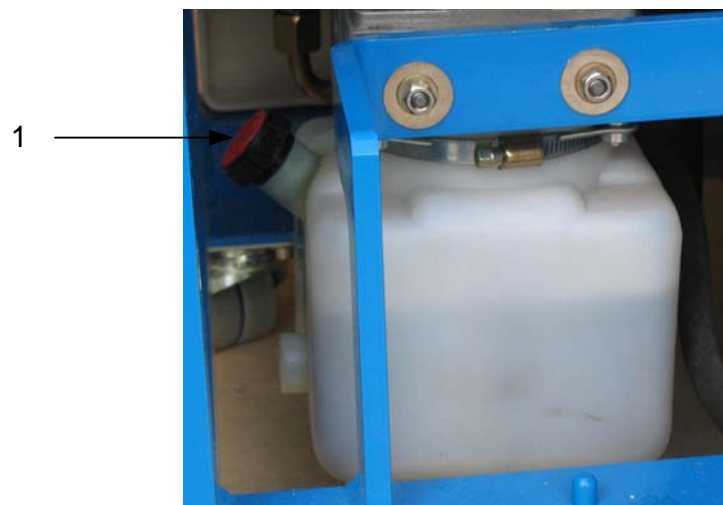


**Figure 7 – Control of oil level**

- Engage the left flank on the locating pins located at the bottom of the tank, press the flank against the chassis and tighten it using the two ¼" fastening screws (Pozidriv no 1 screwdriver).
- Engage the right flank on the locating pins located at the bottom of the tank, press the flank against the chassis and tighten it using the two ¼" fastening screws (Pozidriv no 1 screwdriver).

#### 5.4. Control of cooling liquid level

When the cooling liquid level is too low, an alarm is automatically triggered. If so, you must top up the cooling liquid reservoir with cooling liquid (STERIGERMS must be powered up during this operation):

- Unscrew the two ¼" fastening screws (Pozidriv no 1 screwdriver) on the left flank. Tilt the flank down on your side and pull it upwards to remove it.
- Remove the cap (1) of the cooling liquid reservoir.
- Add demineralised water until the alarm stops.
- Replace the cap (1) of the cooling liquid reservoir.

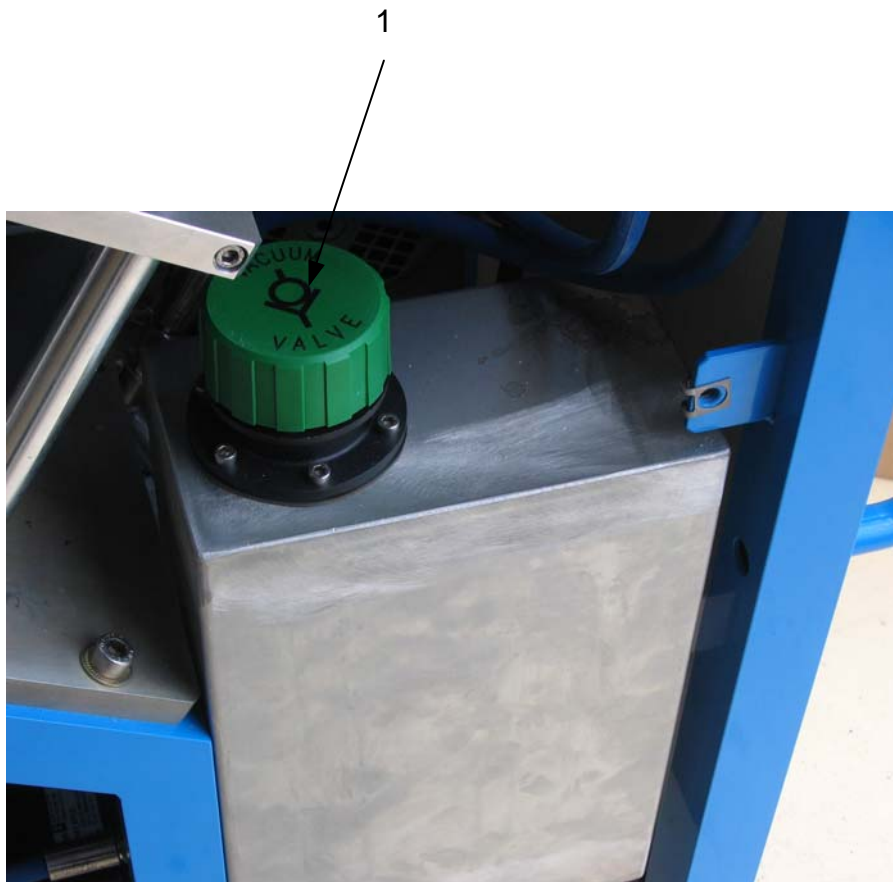

**Figure 8 – Control of cooling liquid level**

- Engage the left flank on the locating pins located at the bottom of the tank, press the flank against the chassis and tighten it using the two ¼" fastening screws (Pozidriv no 1 screwdriver).

### 5.5. Exchange of the cover gasket

- Open the cover of STERIGERMS.
- Remove the gasket (1, Figure 9) using a screwdriver and discard it.

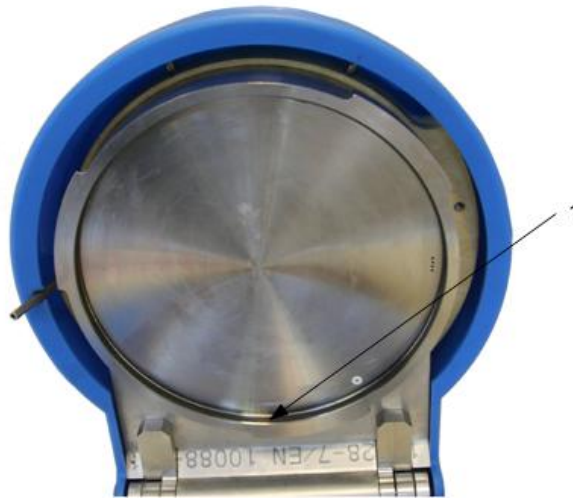


Figure 9 - Exchange of the cover gasket

- Put a new gasket and grease it with **high temperature (180 °C) grease (of type DARINA R2)**.
- Close the cover.


### 5.6. Exchange of the particle filter

- Unscrew the 5 fastening screws (1, Figure 10) (3mm hexagonal key) on the plastic cover of STERIGERMS and remove the plastic cover.



Figure 10 – Removal of the plastic cover

- Unscrew the 2 fastening screws (2 and 7, Figure 11) (3mm hexagonal key) of the filter cap insulating grid (1, Figure 11).
- Open the cover of STERIGERMS and remove the filter cap insulating grid (1, Figure 11).
- Close the cover of STERIGERMS.
- Unscrew the 4 fastening screws (3 and 6, Figure 11) of the filter cap (5, Figure 11) (4mm hexagonal driver).

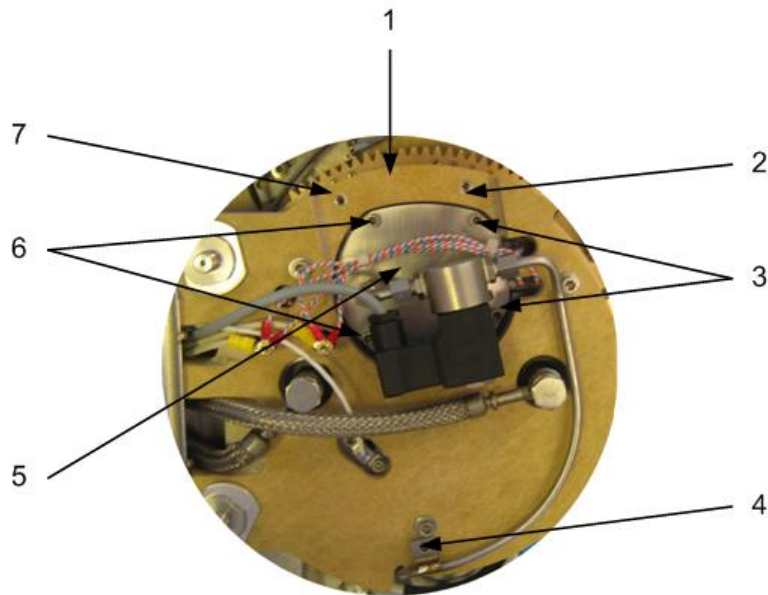

**Figure 11 – Removal of the filter plug**

- Tilt the filter cap (3, Figure 12) to gain access to the particle filter (2, Figure 12) casing.
- Remove and discard the particle filter (2, Figure 12) using a screwdriver.









**Figure 12 – Exchange of the particle filter**

- Clean soiled parts (§ 5.2).
- Check the status of the filter cap gasket (1, Figure 12). You can optionally clean the gasket (§ 5.2) or exchange it if necessary:
  - Remove the filter cap gasket (1, Figure 12) and discard it.
  - Insert a new gasket (1, Figure 12).
- Insert a new particle filter (2, Figure 12) in its casing.
- Put back the filter cap (5, Figure 11) into position and fasten it with its 4 fastening screws (3 and 6, Figure 11) (4mm hexagonal driver).
- Open the cover of STERIGERMS and insert the filter cap insulating grid (1, Figure 11).

- Close the cover of STERIGERMS.
- Screw the 2 fastening screws (2 and 7, Figure 11) (4mm hexagonal key) of the filter cap insulating grid.
- Replace the plastic cover of STERIGERMS and screw the 5 fastening screws (1, Figure 10) (3mm hexagonal key).

## 6. WARNING AND SAFETY SIGNS

	<p><u>Prohibited</u></p> <ul style="list-style-type: none"> <li>- Do not introduce any pressurized containers such as spray bombs.</li> <li>- Do not introduce flammable products.</li> <li>- Do not introduce explosives substances.</li> <li>- Do not perform cycles when the device is empty.</li> <li>- Do not perform prior compaction.</li> <li>- Do not process compressed waste products with a thickness exceeding 10cm, or in other words bags weighing more than 4kg.</li> <li>- Do not introduce any of the products mentioned in the 26 July 1991 circular, volatile toxic substances (including cytostatic drugs) or waste with potential risks of "prion" infection.</li> <li>- Do not introduce waste destined for another treatment line, such as mercury, lead, electric cells, batteries, acids, alkaline compounds radioactive products, etc...</li> <li>- Do not dispose of living animals and/or of recognisable human anatomical parts.</li> <li>- Do not dispose of waste containing plastics for any other purpose than medical or alimentary.</li> <li>- Do not introduce hermetic bags into the device, to guarantee correct operation.</li> </ul> <p><u>Compulsory</u></p> <ul style="list-style-type: none"> <li>- Use the specific Steribag bags developed and validated by the company STERIGERMS and agreed by the "Conseil Supérieur de l'Hygiène Publique de France", indispensable for proper waste processing. Waste products processed without this bag will not be considered as sanitized.</li> <li>- Condition stinging/cutting objects, glass, ceramics and minerals (Sand, gravels) in boxes to that effect. Boxes and other types of waste must be introduced directly into the Steribag bags. Prosthesis and teeth are considered to be cutting/stinging objects.</li> </ul>
	<p>For the compressed waste product to conform to regulations regarding household waste disposal, liquid products cannot exceed 30% of the total weight of a load of processed waste, so as to keep the waste encapsulated in a sufficiently rigid package.</p>
	<p>If waste disinfection is validated by the device (the message "waste processed successfully" appears on the screen), you can dispose of the compressed waste product via the same treatment line as any household waste. If not, the compressed waste product is not disinfected and must imperatively be disposed of via a specialized treatment line for IMW (Infectious Medical Waste) or otherwise reprocessed.</p>
	<p>If there are remaining liquids in the bin after validation by the device of waste disinfection, you can empty the bin from its contents and dispose of it as of any standard liquid household waste product.</p>

	<p>It is forbidden to introduce flammable products.</p>
	<p>It is forbidden to introduce pressurized products such as spray bombs (full or empty).</p>
	<p>Use of the Steribag bag is compulsory. Stinging and cutting objects and glass must be conditioned in a box conform to the norms and intended for this purpose (yellow boxes) before their disposal into the Steribag bag.</p>



## 7. HUMAN – MACHINE DIALOG

STERISOFT management software allows:

- Management of STERIGERMS,
- HMI (Human Machine Interface).

This section introduces available exploitation windows and presents the functions of touch screen keys and related information zones. Unavailable keys are displayed in grey.



TOUCH SCREEN KEYS	FUNCTIONS
Quit	Close the application and shut down the control console with a confirmation window (Window 11)
Enter	Go to the ' <b>Select a user And type your password</b> ' window

**Window 1 - Welcome Please press the “Enter” button**

Stérisoft 12L V4.031c

## Select a user And type your password

User 1  
User 2  
User 3  
User 4  
Administrator

1 2 3  
4 5 6  
7 8 9  
Del. 0 Reset

Validate  
  
Exit

Quit
Enter
23/05/2007 17:40:38

TOUCH SCREEN KEYS	FUNCTIONS
User 1/2/3/4	Select the required user
Validate	Acknowledge and validate the typed password and open the ' <b>Correct password</b> ' or the ' <b>Wrong password</b> ' window. This key is only available when the 4 digits of the password have been typed.
Keyboard	Set of keys allowing to type the password
Exit	Go back to the previous window
Quit	Unavailable
Enter	Unavailable
Del.	Delete the last digit
Reset	Delete all the typed digits

Window 2 - Select a user And type your password

Stérisoft 12L V4.031c

## Correct password

User 1

User 2

User 3

User 4

Administrator

1 2 3

4 5 6

7 8 9

Del. 0 Reset

XXXX

Validate

**Correct password**

Exit

Quit

Enter

23/05/2007  
17:44:36

TOUCH SCREEN KEYS	FUNCTIONS
User 1/2/3/4	Unavailable
Administrator	Unavailable
Validate	Unavailable
Keyboard	Unavailable
Exit	Unavailable
Quit	Close the application and shut down the control console with a confirmation window (Window 11)
Enter	Go to the next window

**Window 3 – Correct password**

StériSoft 12L V4.031c

## Wrong password Please wait and type it again

User 1	1	2	3	XXXX
User 2	4	5	6	Validate
User 3	7	8	9	Wrong password
User 4	Del.	0	Reset	
Administrator				Exit


Quit
Enter

23/05/2007  
17:41:49

Display of this screen is temporary and all the keys are unavailable.


After about 4 seconds, the **'Select a user And type your password'** window appears again, allowing to type the password again.

**Window 4 – Wrong password Please wait and type it again**

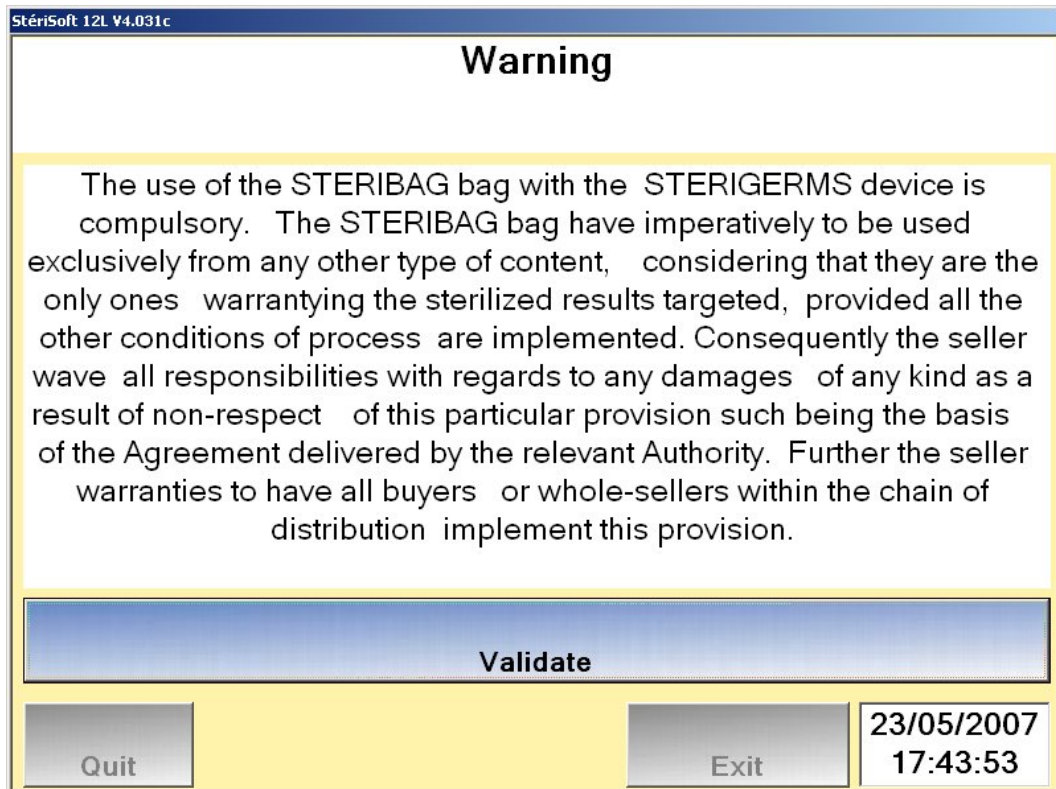


TOUCH SCREEN KEYS	FUNCTIONS
Validate	Go to the next window
Quit	Unavailable
Exit	Unavailable

**Window 5 - Warning (first window)**

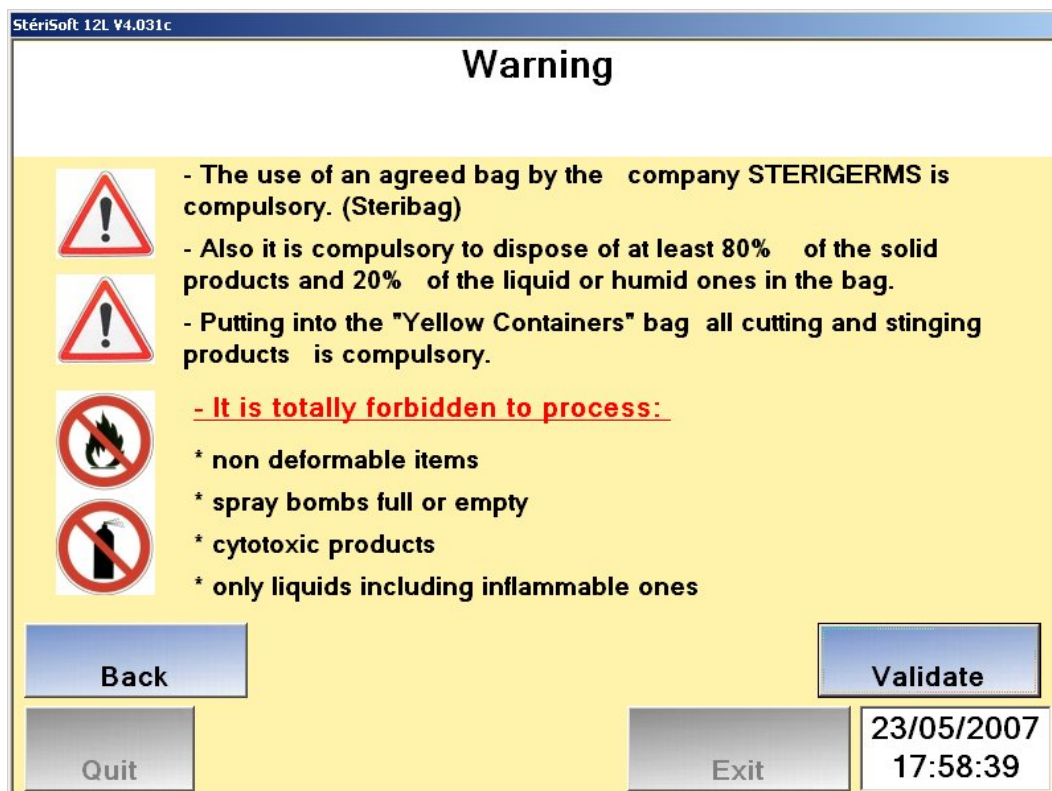

TOUCH SCREEN KEYS	FUNCTIONS
Back	Go back to the <b>'Select an operation'</b> window
Validate	Go to the next window
Quit	Unavailable
Exit	Unavailable


**Window 6 – Warning (second window)**

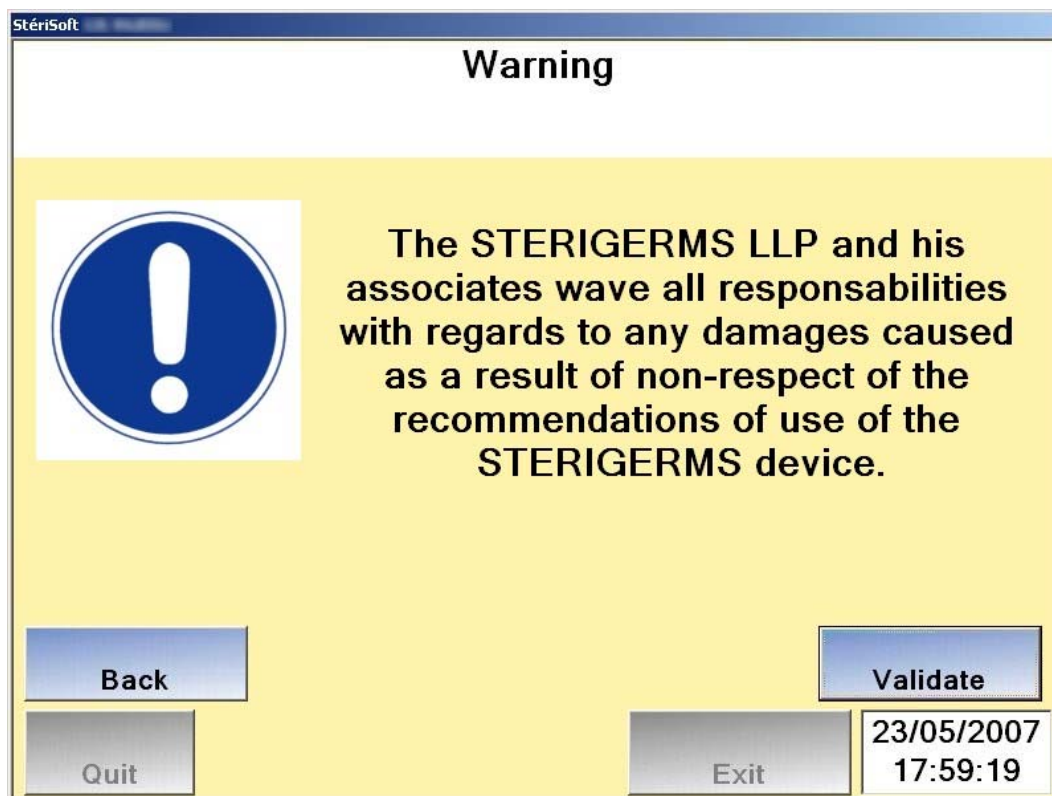

TOUCH SCREEN KEYS	FUNCTIONS
Back	Go back to the second warning window
Validate	Go to the cycle window
Quit	Unavailable
Exit	Unavailable

**Window 7 - Warning (third window)**




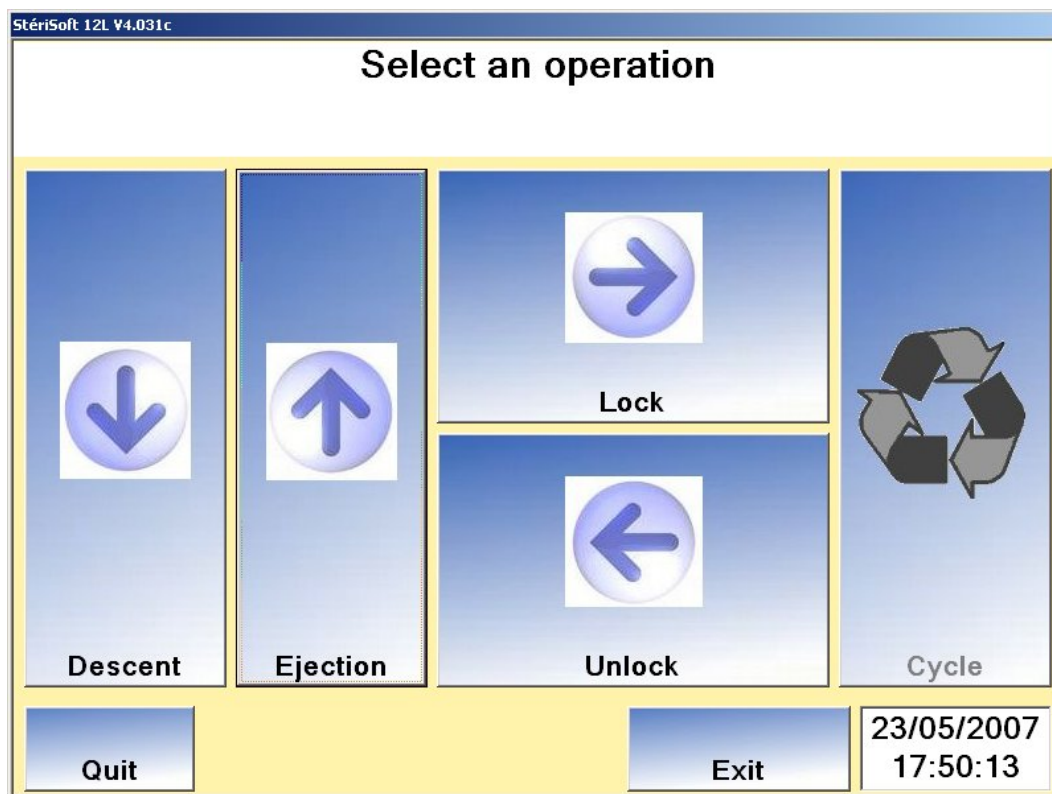

TOUCH SCREEN KEYS	FUNCTIONS
Descent	Control to move the piston downwards
Ejection	Control to move the piston upwards
Lock	Control to lock the cover
Unlock	Control to unlock the cover
Cycle	Go to the cycle window
Quit	Close the application and shut down the control console with a confirmation window (Window 11)
Exit	Go back to the welcome window

#### Window 8 – Select an operation



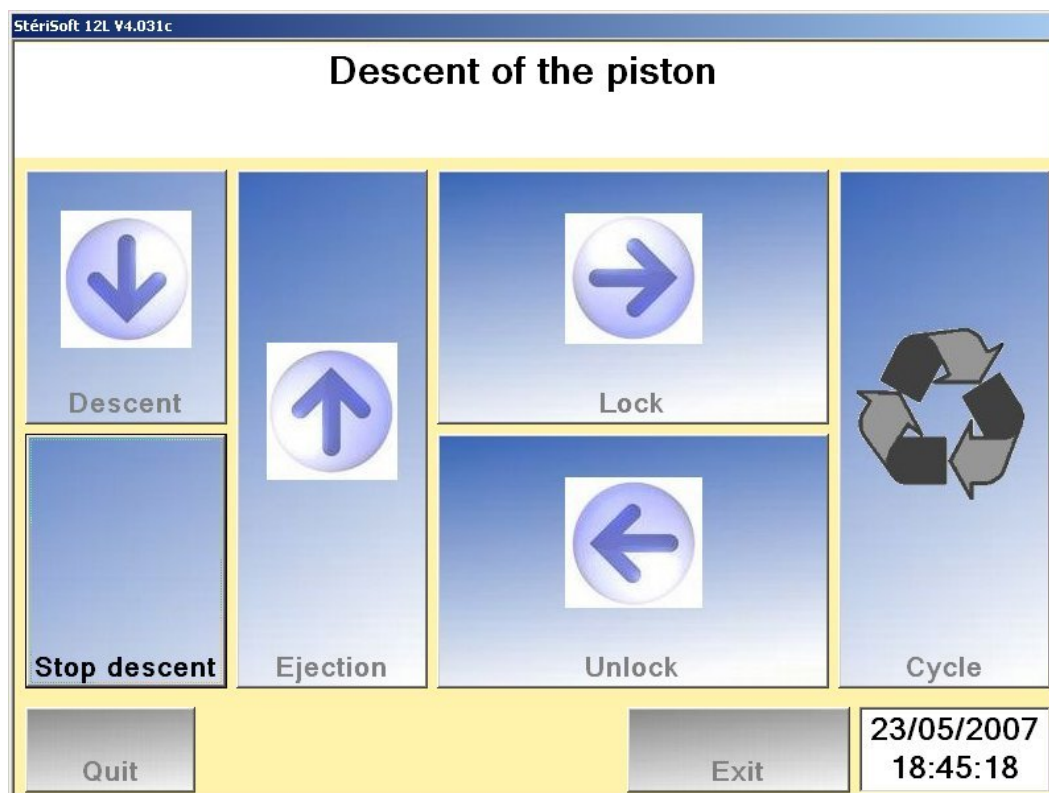


TOUCH SCREEN KEYS	FUNCTIONS
Stop descent	Control to stop the descent of the piston
Ejection	Unavailable
Lock	Unavailable
Unlock	Unavailable
Cycle	Unavailable
Quit	Unavailable
Exit	Unavailable

**Window 9 - Descent of the piston**



TOUCH SCREEN KEYS	FUNCTIONS
Descent	Unavailable
Stop ejection	Control to stop ejection of the compressed waste product
Lock	Unavailable
Unlock	Unavailable
Cycle	Unavailable
Quit	Unavailable
Exit	Unavailable

**Window 10 - Ejection of the compressed waste product**

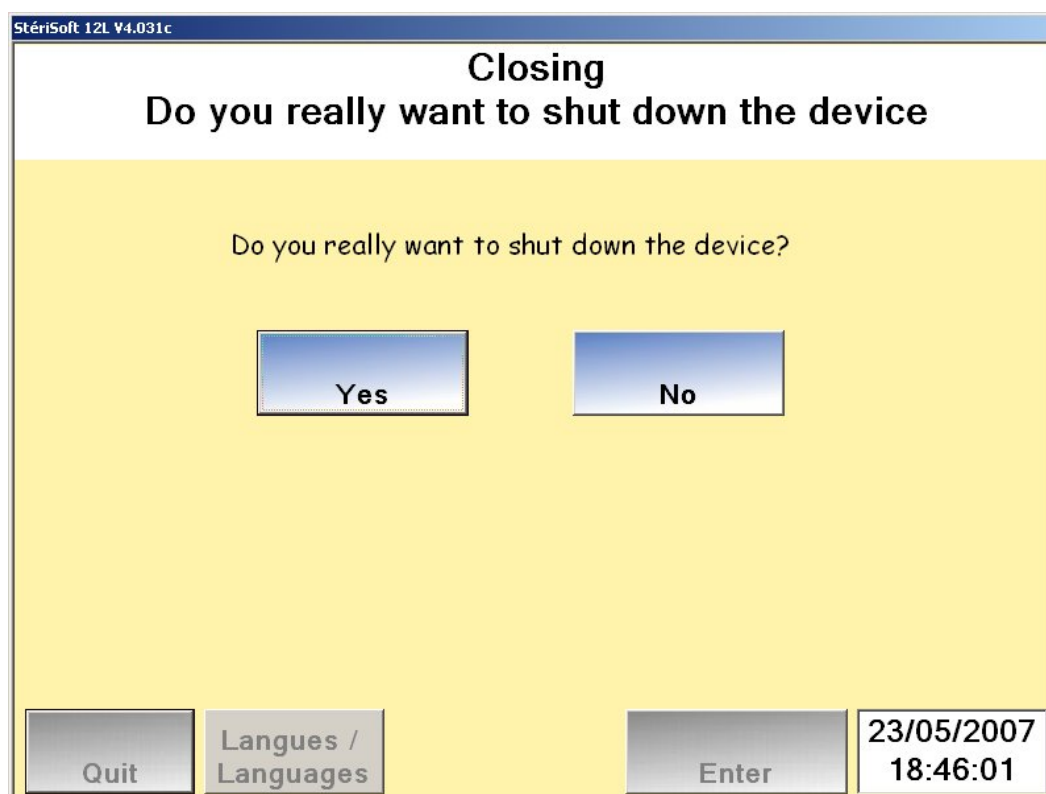

TOUCH SCREEN KEYS	FUNCTIONS
Yes	Close the application and shut down the control console
No	Go back to the previous window
Quit	Unavailable
Enter	Unavailable

**Window 11 - Closing Do you really want to shut down the device**

StériSoft 12L V4.031c

## Press Start to launch the cycle

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


Compression

Heating

Threshold

Cooling

End of cycle



Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product

Start

Back

**23/05/2007**

**18:03:20**

TOUCH SCREEN KEYS	FUNCTIONS
Start	Launch the cycle
Back	Go back to the warning window

INFORMATION ZONES	SIGNIFICATION
Information	Unavailable
Cycle	Unavailable
Compression	Unavailable
Heating	Unavailable
Threshold	Unavailable
Cooling	Unavailable
End of cycle	Unavailable
Current	Electrical consumption of STERIGERMS in sleeping mode
Tension	Tension value of the electrical network
Thick. waste product	Position of the piston (340mm corresponds to low position)

**Window 12 – Press Start to Launch the cycle**


Stérisoft 12L V4.031c

## Press Confirm to launch the cycle

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product



**Confirm**

**Back**

**23/05/2007**  
**18:04:48**

TOUCH SCREEN KEYS	FUNCTIONS
Confirm	Confirm launch of the cycle
Back	Go back to the previous window

INFORMATION ZONES	SIGNIFICATION
Information	Unavailable
Cycle	Unavailable
Compression	Unavailable
Heating	Unavailable
Threshold	Unavailable
Cooling	Unavailable
End of cycle	Unavailable
Current	Electrical consumption of STERIGERMS in sleeping mode
Tension	Tension value of the electrical network
Thick. waste product	Position of the piston (340mm corresponds to low position)

**Window 13 - Press Confirm to launch the cycle**

StériSoft 12L V4.031c

## Starting compression

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


**Compression**

Heating

Threshold

Cooling

End of cycle




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product



23/05/2007

18:30:21


INFORMATION ZONES	SIGNIFICATION
Information	Display of detected malfunction (§4)
Cycle	Cycle number
Compression	Display of compression phase
Heating	Unavailable
Threshold	Unavailable
Cooling	Unavailable
End of cycle	Unavailable
Current	Electrical consumption of STERIGERMS in sleeping mode
Tension	Tension value of the electrical network
Thick. waste product	Position of the piston (340mm corresponds to low position)

**Window 14 – Start compression**

Stérisoft 12L V4.031c

## Heating

**Information**  
 Liquid detected  
 Waste product too thick  
 Power outage  
 Low tension  
 Pressure too high  
 Cold liquid leakage



Cycle: 1


Compression

Heating  
23/05/2007 18:33:34

Threshold

Cooling

End of cycle

  
 Version 12L V4.031c

Current (A) 1  
 Tension (V) 230  
 Thick. waste product 53

23/05/2007  
18:33:37

INFORMATION ZONES	SIGNIFICATION
Information	Display of detected malfunction (§4)
Cycle	Cycle number
Compression	Display of completed compression phase
Heating	Display of heating phase parameters
Threshold	Unavailable
Cooling	Unavailable
End of cycle	Unavailable
Current	Electrical consumption of STERIGERMS during heating phase
Tension	Tension value of the electrical network
Thick. waste product	Position of the piston (340mm corresponds to low position)

**Window 15 - Heating**

StériSoft 12L V4.031c

## Threshold reached

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


**Compression**

**Heating**  
23/05/2007 18:35:55

**Threshold remain 12**  
23/05/2007 18:35:55

**Cooling**

**End of cycle**




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product



23/05/2007  
18:36:01

INFORMATION ZONES	SIGNIFICATION
Information	Display of detected malfunction (§4)
Cycle	Cycle number
Compression	Display of completed compression phase
Heating	Display of heating phase parameters
Threshold	Display of processing phase parameters
Cooling	Unavailable
End of cycle	Unavailable
Current	Electrical consumption of STERIGERMS during processing phase
Tension	Tension value of the electrical network
Thick. waste product	Thickness of the compressed waste product

**Window 16 – Start threshold**



Stérisoft 12L V4.031c

## Cooling

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


**Compression**

**Heating**  
23/05/2007 18:37:30

**Threshold**  
23/05/2007 18:37:30

**Cooling**  
23/05/2007 18:37:30

**End of cycle**




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product



23/05/2007  
18:37:37

INFORMATION ZONES	SIGNIFICATION
Information	Display of detected malfunction (§4)
Cycle	Cycle number
Compression	Display of completed compression phase
Heating	Display of heating phase parameters
Threshold	Display of processing phase parameters
Cooling	Display of cooling phase parameters
End of cycle	Unavailable
Courant	Electrical consumption of STERIGERMS during cooling phase
Tension	Tension value of the electrical network
Thick. waste product	Thickness of the compressed waste product

**Window 17 - Cooling**

StériSoft 12L V4.031c

## End of Cycle

**Information**

- Liquid detected
- Waste product too thick
- Power outage
- Low tension
- Pressure too high
- Cold liquid leakage

Cycle:


**Compression**

**Heating**  
23/05/2007 18:41:03

**Threshold**  
23/05/2007 18:41:03

**Cooling**  
23/05/2007 18:41:03

**End of cycle**  
23/05/2007 18:41:03




Version 12L V4.031c

Current (A)

Tension (V)

Thick. waste product



Waste has been processed

Opening device

23/05/2007  
18:41:08

TOUCH SCREEN KEYS	FUNCTIONS
Open the device	Go back to the ' <b>Select an operation</b> ' window

INFORMATION ZONES	SIGNIFICATION
Information	Display of detected malfunction (§4)
Cycle	Cycle number
Compression	Display of completed compression phase
Heating	Display of heating phase parameters
Threshold	Display of processing phase parameters
Cooling	Display of cooling phase parameters
End of cycle	Display of end of cycle phase parameters
Waste has been processed	Display of processing result after cycle completion
Current	Electrical consumption of STERIGERMS after cycle completion
Tension	Tension value of the electrical network
Thick. waste product	Thickness of the compressed waste product

**Window 18 – End of cycle**

StériSoft

## End of Cycle





**Cleaning of the machine is a compulsory rule after every use**

Gasket

Cover

Airholes

Recovery tray

Cover	Clean with a wet sponge the 4 elements
Piston	
Gasket	
Airholes	
Recovery tray	Empty regularly the recovery tray

**Non compliance with this rule entails the guarantee cancellation**

OK

23/05/2007  
18:41:08

TOUCH SCREEN KEYS	FUNCTIONS
OK	Go back to the 'Select an operation' window

### Window 19 – Cleaning of the machine



Window 20 - Machine error Only administrator can enter

## **8. WARRANTY**

The STERIGERMS contractual warranty shall not apply in case of non-observance of any of the instructions as set out in the User's Manual.

The device complies with the following European directives:

89/392/CEE dated 06/14/89 modified

73/23/CEE dated 02/19/73 modified

89/336/CEE dated 05/03/89 modified

